



# ONEFILL™ Carrier-Based Obturators

**Easy-to-use plastic obturators feature a new level of precision, performance, and radiopacity!**

New **OneFill™** can be used to obturate any constant or variable 04, 06, 08, and 10 tapered canal, size 20 to 50, shaped by the SAF Endo System, NiTi rotary systems, and hand-filing techniques.

- Heated alpha-phase gutta percha provides excellent 3-D fill
- Precise coating for improved accuracy and easier cleanup
- Popular ISO sizes in 25mm lengths for system versatility
- Extremely radiopaque and easy to use
- Compatible for use in leading oven brands
- Calibration rings designed for easy working length determination
- New Size Verifiers match obturators and confirm precise sizing



## ONEFILL™ Carrier-Based Obturators Ordering Information



<b>OneFill Obturator Refill Packages .....</b>	<b>20 per Pkg.</b>	<b>\$139.99</b>
		5 Pkg. @ \$129.99
Specify Size:	20 (704-0036)    30 (704-0038)    40 (704-0040)    50 (704-0042)	
	25 (704-0037)    35 (704-0039)    45 (704-0041)	



<b>OneFill Obturator Refill Packages .....</b>	<b>6 per Pkg.</b>	<b>\$49.99</b>
		5 Pkg. @ \$43.99
Specify Size:	20 (704-0034)    30 (704-0029)    40 (704-0031)    50 (704-0033)	
	25 (704-0035)    35 (704-0030)    45 (704-0032)	



<b>Size Verifier Packages .....</b>	<b>6 per Pkg.</b>	<b>\$19.99</b>
Specify Size:	20 (704-0043)    30 (704-0045)    40 (704-0047)    50 (704-0049)	
	25 (704-0044)    35 (704-0046)    45 (704-0048)	

**Find out more about this new product • [www.hsendo.com](http://www.hsendo.com)**

Exclusively distributed by  
**HENRY SCHEIN**<sup>®</sup>  
 DENTAL

To Order: **1-800-372-4346** 8am-9pm (et)  
 To Fax: **1-800-732-7023** 24 Hours  
 Pricing expires December 30, 2011

TM-The OneFill brand name is owned by US Endo and Exclusively used by Henry Schein, Inc.

Easy-to-use plastic obturators feature a new level of precision, performance, and radiopacity!


**OneFill™ carrier-based obturators are easy to use. Simply choose the OneFill obturator ISO size that matches the largest ISO file size taken to length in the canal.**

*For example, if you finish the canal with a Size 20 file, with an 04, 06, 08, or 10 taper, then fill with the Size 20 OneFill obturator. Use the OneFill tip size matching the tip size of your largest file. Additionally, OneFill can be used for any tapered canal created by any instrument, including the SAF system, NiTi rotary systems, and hand filing.*




### Step-by-Step Instructions

- Place the OneFill obturator in a 5.25% sodium hypochlorite solution for one minute, then rinse with 70% alcohol and allow carrier to dry thoroughly before use.

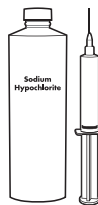

- Insert the OneFill obturator into the oven and start the heating cycle.

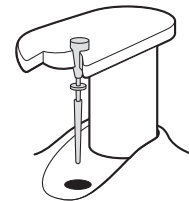


- Dry the canal and apply a thin layer of sealer.



- Remove the obturator when the oven beeps.




- Set the stopper to the working length (or use the calibration rings).

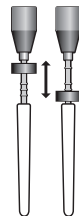


- Immediately place the obturator smoothly in the canal to the working length.



- OneFill can be used with the same heat setting as your current oven. Set the oven to the corresponding ISO size.

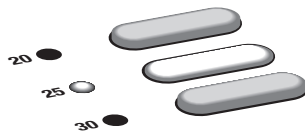

- Cut the obturator at the orifice level.

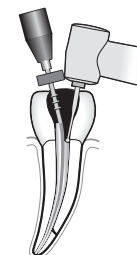

Find out more about this new product • [www.hsendo.com](http://www.hsendo.com)

Exclusively distributed by

**HENRY SCHEIN®**  
DENTAL

To Order: **1-800-372-4346** 8am-9pm (et)

To Fax: **1-800-732-7023** 24 Hours

Pricing expires December 30, 2011